

Your result is displayed in units called ACEs.

After displaying your measurement, Biosense will automatically go to sleep.



0-4	No or Low Ketones
5-14	Moderate Ketosis
15-30	Advanced Ketosis
31-39	Deep Ketosis
40+/HI	High Ketosis

The color scale indicates where your measurement lies on the spectrum of ketosis.

Learn more about the ACEs unit:
mybiosense.com/understanding-the-aces-unit

Tips for Success:

Store Biosense in the zippered pouch when not in use. Be sure to keep Biosense away from solvents such as:

- hand sanitizer
- cleaning solution
- nail polish remover
- perfume

Wait 60 minutes to use Biosense after eating, smoking, vaping, drinking flavored beverages, or using products such as:

- toothpaste
- mouthwash
- cough drops
- breath mints
- chewing gum
- lip balm
- artificial sweeteners


After consuming alcohol, wait until the next day to use Biosense.


Frequently Asked Questions: mybiosense.com/faq

Questions? Email us at support@mybiosense.com

Readout, Inc.
2940 Locust Street
St. Louis, MO 63103, USA

Stay up to date :

 @mybiosense

 @mybiosense

 mybiosense.com





Breath Ketone Monitor



Quick Start Guide


The first time you use Biosense, and later if Biosense loses power, you will be prompted to pair to the mobile app to set the device clock.

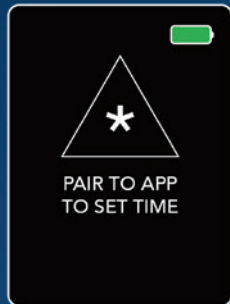

If you choose to use Biosense without the app, press  to skip pairing. If you skip pairing, Biosense will prompt you to set the device clock manually.


*For directions to pair to the app, visit:
mybiosense.com/pairing*

How to take a measurement:


To wake up, press and hold the  button until the logo appears, then release the button.




Press START  to begin a measurement.

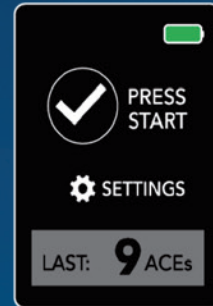


After a normal inhale, place your lips around the mouthpiece and exhale with steady medium force.

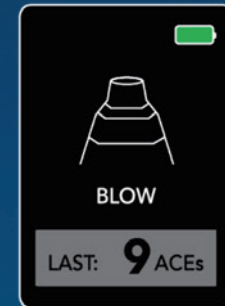

Biosense will take up to 300 seconds to warm up.

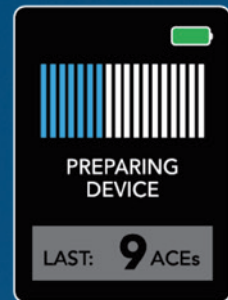

Biosense will make final preparations.

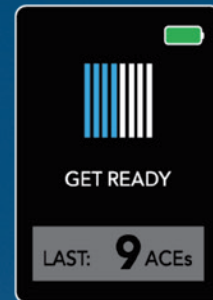

Do not stop blowing until after you feel Biosense stop vibrating.

Biosense will take about 15 seconds to analyze your breath sample.

